


The highest performance of sample separation achieved with the use of Universal Column Premium (30um Silica Gel) and Inject Column.




Synergistic effect of using Universal Column Premium and Inject Column in pairs

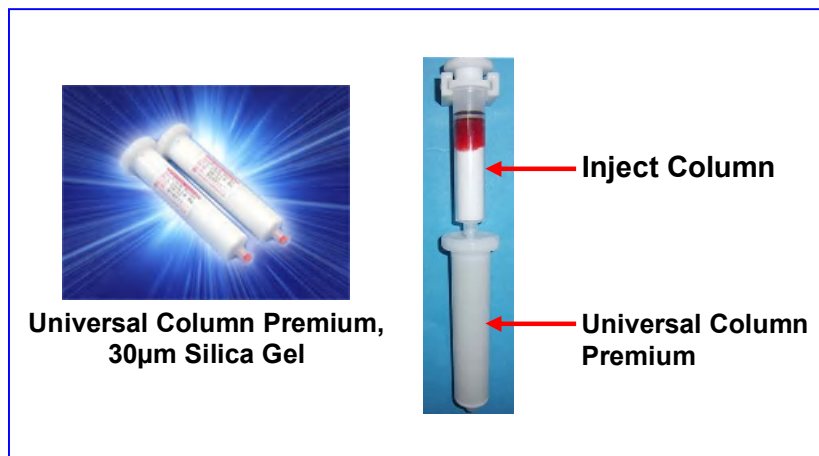

Hexane/Ethyl acetate = 70/30
Butyl-p-hydroxybenzoate
Methyl-p-hydroxybenzoate

Premium Universal Column L-size and Inject Column M-size

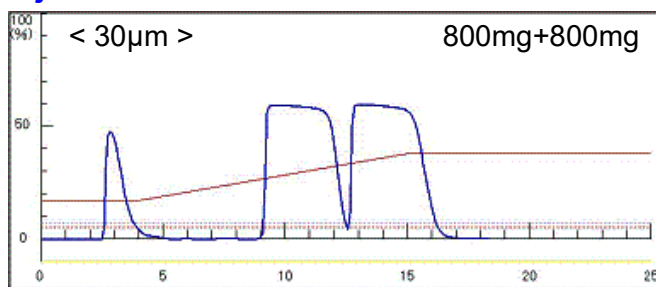


Premium Universal Column L-size
(30 μ m silica gel, 40g, I.D.2.6 x 12.4cm packed length, O.D.3.0 x 16.5cm overall length)




Sample load:
2 times bigger

Premium Universal Column L-size, Direct Sample Loading

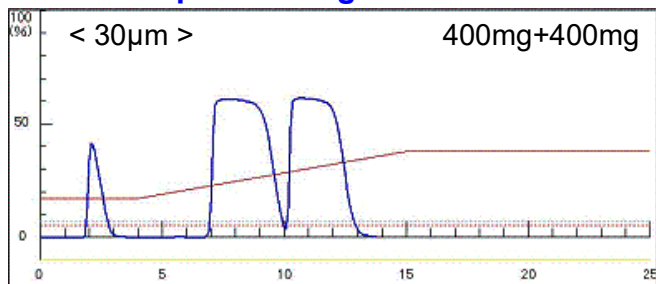


Premium Universal Column L-size
(30 μ m silica gel, 40g, I.D.2.6 x 12.4cm packed length, O.D.3.0 x 16.5cm overall length)

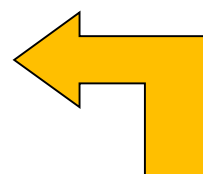


Sample load:
4 times bigger

Standard Universal Column L-size, Direct Loading

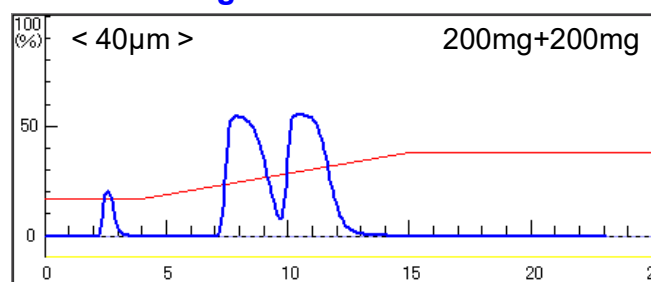

Standard Universal Column L-size
(30 μ m silica gel, 40g, I.D.2.6x12.4cm packed length, O.D.3.0 x 16.5cm overall length)

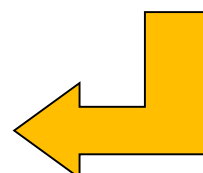

Sample load:
2 times bigger